

E8

EN

Efficient and reliable,
sterilisation for everyone



Euronda[®]
Pro System

E8

Smart sterilisation is a question of balance



Reliability

Control systems and vanguard functions make sure that processes and results are consistent over time a Euronda Pro System hallmark.



Performance

Powerful, fast and smart: it keeps consumption down and adapts cycle length to load.



Design

Elegant, compact and modern: E8 an harmonious combination of fine aesthetics and total functionality.



Ergonomics

Easy to install and even easier to use. E8 has been designed to save space and make healthcare professionals' lives simpler.



Innovation

All Euronda Pro System's top technology concentrated in an autoclave that suits any practice.



Traceability

E8's traceability system allows for saving data on SD Card and for the integration two different kinds of printer, according to needs.

Straightforward, high-performing and affordable: E8 brings together all Euronda's top technology to strike a perfect balance between innovation & ergonomics, power & consumption and customisation & usability. The end result is a state-of-the-art autoclave which perfectly meets the requirements of all practices. E8 boosts efficiency and makes your job easier, whilst delivering total reliability and longevity. E8 is the most affordable way to reap the advantages of all the Pro System sterilisation solutions.

Reliability

Efficiency and safety that stand the test of time

Constant safety and steadfast efficiency - this is what E8 is all about. What you get is a smart sturdy autoclave of the utmost flexibility developed to adapt to the various requirements of any practice and to outperform. E8 integrates and harmonises the best solutions and functions of all Euronda Pro System devices. The various components, control systems and software strive together actively to give you absolute reliability.

Euronda

Pro System

E8





Water quality sensor

The built-in conductivity meter monitors water quality and, if it is not satisfactory, E8 informs the operator with special messages on the display.



Dirt Control System

It is located in the hydraulic circuit and it makes the dirt easily flow in the used-water tank, so that the autoclave can work without problems.



Air-water separator

The air-water separator device, designed and built into the tank unit, is a real “life saver” for the vacuum pump. The pumps used on machines of this level are essential to achieve proper sterilisation, but are very delicate as they are in constant contact with flows containing high temperature steam, air, dirt and chemical waste from lubricating oils found on instruments being sterilised (wands, turbines,...). These streams are channelled into the separator device which separates the air from the other components before they enter the vacuum pump. This phase allows the pump to always operate in ideal working conditions, preserving the efficiency of the device.



Support software

Euronda has developed a number of software applications to support technicians during the installation, diagnosis and maintenance processes. This helps to guarantee peak performance from E8 at all times while cutting assistance times and costs.



The instant steam generator and enhanced drying cycle ensure a rapid and superlative sterilisation process which keeps down water (300 ml of water per cycle) and energy (power usage of 2300 W at 10 A) consumption by adapting to the load.



Performance

The ideal balance between power and consumption



Efficient

The steam generator instantly produces the steam for the sterilisation chamber, optimizing the water consumption in accordance with the load, making further cuts to the sterilisation time and actually improving the quality of the steam generated.



Customization

First and foremost, E8 can run five different type B cycles (121°C, 134°C, 134°C Prion, 134°C Rapid and 134°C Rapid Prion) in full compliance with the EN 13060 standard, as well as three test cycles (Vacuum, Helix, and Bowie and Dick). In addition, thanks to the password-activated Personal functions, it is also possible to run Light cycles (type N cycles - N134 and N121) and Light & Stock cycles (type S cycles - S134 and S121).



Eco-friendly

Fast and efficient at the same time. The Euronda Pro System technology optimizes the sterilisation settings to guarantee savings in terms of water (300 ml of water per cycle), energy (power usage of 2300 W at 10 A) and time.

Design

A perfect blend of aesthetics and functionality

With its elegant contemporary look, E8 underpins and reinforces the new style launched by Euronda autoclaves. The Soft Touch colour display is protected by a polycarbonate screen and enhanced with a distinctive front panel in light grey.



Elegant

Clean modern lines for a practical reliable autoclave endowed with all the signature details of Euronda quality.



Functional

The exclusive Soft Touch display comes with a protective polycarbonate screen designed for easy cleaning.



Compact

Its linear design, compact size and adjustable supports mean it easily fits into small spaces or on the work top.



Innovation

Top technology at your fingertips

E8 encapsulates the best of Euronda's Pro System cutting-edge technology and brings it within everyone's reach. Programming and customisation functions, an extremely user-friendly interface and incredibly innovative systems to keep the autoclave's performance standards consistent over time are all technological features which in E8 are synonymous with simplicity.



The smart autoclave

Readily understandable and easy to use, E8 is a prime example of how simple and intuitive technology can be.



Soft Touch display

The Soft Touch colour display and brand-new clear fuss-free icons enable healthcare professionals to work intuitively and fast. The user interface has been specially designed to speed up and optimise routine work in the practice.



E-Timer

Make the most of your time! Thanks to the innovative E-Timer, you can set tests and sterilisation cycles to start at whatever time you like, allowing you to save time and start the working day at full speed.

Ergonomics

E8 makes everything easier



Devised to fit into your space

Endowed with a new rear ventilation system, it fits snugly on the shelves of the sterilisation room.



Energy efficiency at your service

The choice of containing the power used within 2300 W allows the installation on a 10 A line.



Inspection[®] System

The tank compartment access system is extremely quick and simple. The upper cover is made of an innovative engineering plastic. It is held in place by a magnetic system which makes it easy to open at any time and enables even less proficient people to clean it easily and frequently. This guarantees safe sterilisation that remains consistently effective in the long term.



Modular structure

The side panels are made of cutting-edge materials and they are extremely simple to remove. In addition, the components inside the machine are easily accessible and positioned in a logical manner. This helps to make maintenance easier, faster and more effective.



Tank compartment access system

Traceability

Keep the sterilisation process constantly under control

Traceability is a hallmark of all Euronda Pro System autoclaves, the first to have an expandable traceability system. E8 is no exception to the rule. All sterilisation cycles are saved to an SD Card for transferral onto a computer. With E8, two different kinds of printers can be added according to your requirements, so you can save time without compromising safety standards.



Print Set 1
Built-in thermal
paper printer



Print Set 2
Built-in adhesive
label printer

Euronda
Pro System

E8



Expandable connectivity

A storage system on a standard integrated SD Card so that all data can be easily transferred from the autoclave to the computer, in compliance with EN 13060 standards.



Validation and load release

This function allows users to keep a digital record of the compliance of the sterilisation process in both the loading phase and the unloading phase. The operator ensures that the pouches are dry and checks their integrity and the state of the instruments, while also making sure that the sterilisation parameters are satisfactory.



Optional printers

E8's expandable traceability system caters to a range of needs, with a selection of two different optional printers. E8 can satisfy the traceability requirements of any practice and guarantee compliance with the legally established safety standards.



Starter kit included



Standard tray holder

Square profile steel tray holder with a capacity of five standard aluminium trays. If rotated by 90°, it can hold three sterilisation cassettes.

Dimensions (wxhxd):

186 x 22 x 379 mm

Five trays included



- Drainage pipes
- Sponge
- Tray removal handle
- Door adjustment handle
- Funnel

Sterilisation cycles

	B 134	B 121	
	Solid, porous and hollow, wrapped and unwrapped	Solid, porous and hollow, wrapped and unwrapped. Particularly indicated for plastics, turbines and handpieces	
Temperature	134° C	121° C	
Pressure	2.05 bar	1.05 bar	
Sterilisation time	4'	20'	
Drying time	15'	15'	
Total time with standard loads*	29'	44'	
Maximum load solid/porous	6 Kg / 2 Kg	6 Kg / 2 Kg	

Personal functions

(activation on demand with a password)

	Light & stock cycles**		
	S 134	S 121	
	Sterilisation of non-bagged or bagged solid products	Sterilisation of non-bagged and plastic or bagged solid products	
Temperature	134° C	121° C	
Pressure	2.05 bar	1.05 bar	
Sterilisation time	4'	20'	
Drying time	15'	15'	
Total time with standard loads*	18'	32'	
Maximum load solid	6 Kg	6 Kg	

* Times are intended without drying. The times have a tolerance of +-5%.

** All of the cycles are available unless it is not possible due to specific national legislation.

	B 134 prion	B 134 - rapid	B 134 prion - rapid
	Prion cycle for wrapped and unwrapped solid, porous and hollow	Rapid cycle for wrapped and unwrapped solid, porous and hollow For 1 tray instruments	Rapid prion cycle for wrapped and unwrapped solid, porous and hollow For 1 tray instruments
	134° C	134° C	134° C
	2.05 bar	2.05 bar	2.05 bar
	18'	3'5"	18'
	15'	4'	4'
	43'	27'	42'
	6 Kg / 2 Kg	0.6 Kg / 0.2 Kg	0.6 Kg / 0.2 Kg

Light cycles**	
N 134	N 121
Unwrapped solid instruments	Sterilisation of non-bagged and plastic solid products
134° C	121° C
2.05 bar	1.05 bar
4'	20'
4'	4'
24'	33'
6 Kg	6 Kg

Technical data

	E8 24L
Overall dimensions	460 x 455 x 610 mm (wxhxd)
Weight (machine empty with basket)	48 Kg
Power supply	230 V 50/60 Hz
Absorbed current	2300 W 10 A
Operating temperature	+5 - +40°C
Operating altitude	≤ 2000 m a.m.s.l.*
Minimum depth of the support surface	500 mm
Minimum width of the support surface	460 mm
Minimum space above the tank	0 mm leaving the loading cap free
Noise level	50 dB
Weight per support area (tank full and chamber with max. load)	48 + 9 = 57 Kg Load on each foot: 14.3 Kg
Operating range	Indoor use - Maximum relative humidity 85% Maximum line-voltage deviation ±10% Installation category II - Pollution level 2
Average water consumption per cycle	300 ml
Tank capacity (clean water / used water)	5 / 3.5 litres

* For altitudes higher than 2000m, please contact Euronda SPA

Safety Device

- Overpressure protection
- Accidental door opening protections
- Process appraisal system
- Real-time self diagnostic system

Compliance with harmonised standards

EN 13060
EN 14971
EN 61010-1
EN 61010-2-040
EN 61326
EN 13445

Documentation


- User Manual
- Declaration of conformity chamber
- Declaration of conformity autoclave
- Warranty Certificate
- Installation chart
- Servicing handbook

Certification information

- E8 Autoclave compliant with Directive 93/42/EC and addenda.
- Class II b device
- Stainless steel moulded chamber certified to PED Directive 2014/68/EC
- Chamber capacity 24L
- 5 trays o 3 cassettes
- Company certification: ISO 9001, ISO 13485
- CE Type Ref: E8 24L

Technical Specifications

Overall required space (mm)

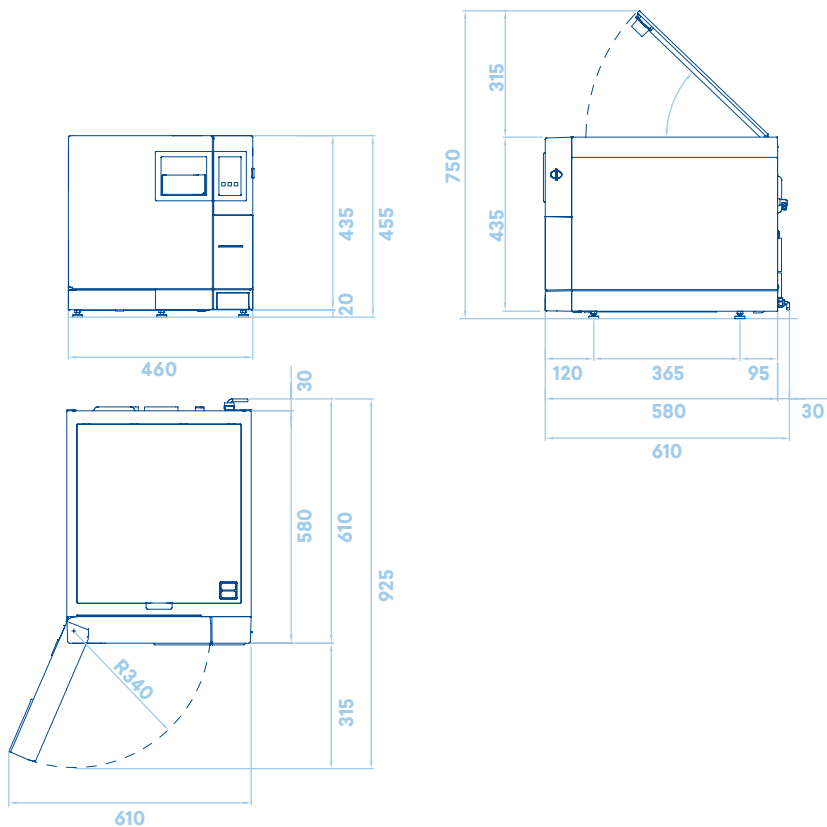

Sterilisation chamber sizes (mm)

